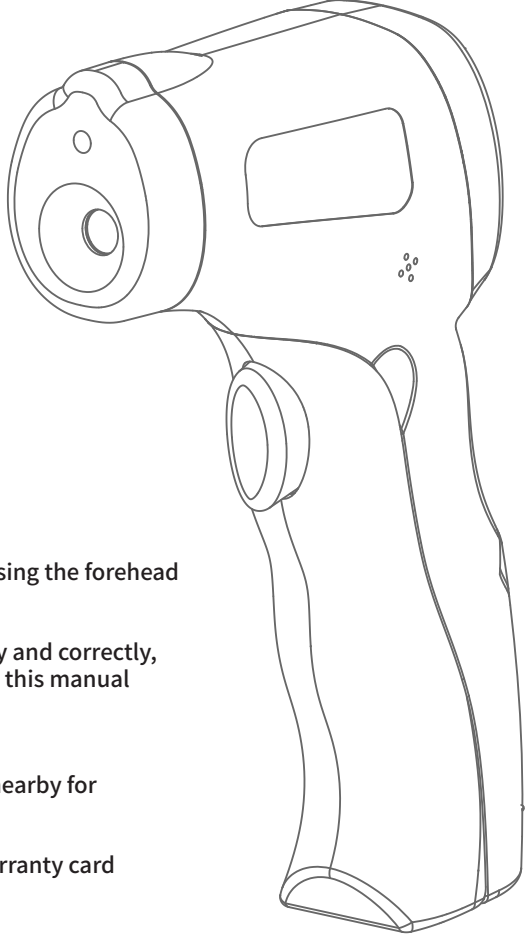


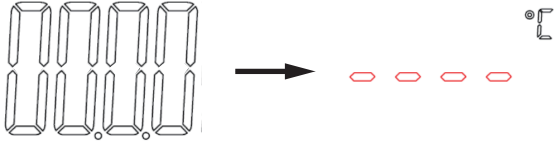
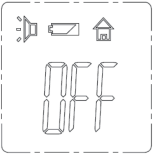





正面

<div><div>MEDICAL INFRARED THERMOMETER</div><div></div><div><p>Thank you very much for purchasing the forehead infrared thermometer.</p><p>In order to use this product safely and correctly, please read and fully understand this manual</p><p>Safety precautions for.</p><p>Please always keep the manual nearby for reference.</p><p>(please keep this manual and warranty card properly and do not lose it.)</p></div></div>	<div><div>Precautions Warning</div><div><p>1. It is very dangerous for patients to judge and treat themselves only by measuring results, so please follow the instructions of doctors.</p><p>● self judgment may lead to deterioration of the condition.</p><p>2. Do not touch or blow the infrared sensor with your hand or mouth, or let the sensor face the steam.</p><p>● when the infrared sensor is damaged or dirty, it may cause abnormal measurement results.</p><p>3. If there is any temperature difference between the storage site and the measuring site, please leave it at room temperature (measuring site) for about 30 minutes before measuring.</p><p>● may cause abnormal measurement results.</p><p>4. Please keep the machine out of the reach of children.</p><p>● if you accidentally eat the battery or battery cover, please contact your doctor immediately.</p><p>5. When measuring the temperature, please do not close to the position of the heating and cooling or the air outlet.</p><p>● to avoid affecting the accuracy of measurement.</p><p>6. When the infrared sensor is dirty, please use cotton ball or soft dry cloth dipped in 95% anhydrous alcohol to gently wipe.</p><p>● if you wipe it with toilet paper or paper towel, it will scratch the infrared sensor and lead to incorrect measurement results.</p><p>7. Mechanical damage such as heavy fall or collision shall be avoided.</p><p>● it is possible that the measurement is inaccurate.</p></div></div>	<div><div><p>1. The distance between the front temperature gun and the forehead must be less than 5cm.</p><p>● to avoid affecting the accuracy of measurement.</p><p>2. Do not drink, eat or exercise before or during measurement, and do not measure when sweating.</p><p>● it is likely to affect the accuracy of measurement.</p><p>3. Do not discard the battery to a dangerous area or discard it at will.</p><p>● the battery may break.</p><p>● pollute the environment.</p><p>4. This product has self-test function.</p><p>5. Cosmetics and skin color, due to different infrared radiation rate, will affect the display temperature accuracy.</p><p>● if the person with black skin needs to manually adjust the temperature compensation value before measurement.</p><p>6. Do not soak in water or direct sunlight.</p><p>7. When measuring the forehead temperature of human body, the product must be operated in the forehead temperature mode (see point 10 of operation steps).</p></div><div><div>Proposal</div><div><p>1. When you tell the doctor the temperature you have taken, please indicate that you have taken it with a forehead infrared thermometer.</p><p>2. Please do not force collision, fall, trample or shake the product.</p><p>3. Do not disassemble, repair or modify this product.</p><p>4. The product is not waterproof. Please take care not to let liquid (alcohol, water drop, hot water, etc.) enter into the product.</p><p>5. The product must be kept clean and stored in a dry place.</p></div></div></div>
---	--	---

反面

<div><div><p>6. If you find any problems, please contact the sales department. You can't repair the product by yourself.</p><p>7. Frontal temperature gun is vulnerable to radiation interference. Please do not use it in the environment with strong electromagnetic interference.</p><p>8. Please dispose the waste and residue at the end of the service life of the product according to the local laws and regulations.</p></div><div><div>Usage method</div><div><p>1. After opening the package, open the battery door from "figure 1", and install two AA batteries No.5 as shown in "Figure 2" (pay attention to the positive and negative poles of the battery)</p><div><div></div><div></div></div><div><div>figure 1</div><div>Figure 2</div></div><p>2. After the battery is installed, press the measurement key once to start the machine. The four "0" numbers in the middle of the LCD display become "----", as shown in Figure 3.</p><div></div><div><div>Figure 3</div></div><p>Press the measurement key once, LCD will display "FB" → "LD" → "temperature", FB = object LD = human body</p></div></div></div>	<div><div><p>4. If there is no operation for 30 seconds, the LCD screen will display off, and it will shut down automatically after 5 seconds.</p><div></div></div><div><p>5. In the current temperature measurement data display mode, press the "△" key to turn on the radium spotlight function.</p></div><div><div><p>6. In the current temperature measurement data display mode, press "M" to enter the "object" temperature measurement or "human body" temperature measurement display mode,</p><div><div></div> = “Human body”<div></div> = “Object”</div></div><div><div></div></div></div></div>	<div><div>Maintenance and safekeeping</div><div><p>After use, please maintain and store this product.</p><p>1. Please use a soft dry cloth to gently wipe the dirt on the body.</p><p>Do not use water or detergent, diluent or volatile oil containing abrasives.</p><p>2. When the infrared sensor is soiled, please use a soft dry cloth or cotton ball dipped in 75% anhydrous alcohol to gently wipe it. If it still cannot be wiped dry</p><p>Please contact the service center when cleaning.</p><p>Do not use toilet paper or tissue to wipe.</p><p>3. Please do not store in the place with direct sunlight, high temperature and humidity, much dust, close to fire and vulnerable to vibration and shock.</p><p>4. When the forehead temperature gun is not used for a long time, take out the battery to avoid battery leakage.</p><p>5. When the battery symbol "" appears on the screen of the forehead temperature gun, it indicates that the battery power is insufficient and needs to be replaced as soon as possible.</p><p>6. This product adopts special chip and sensor, with stable performance and reliable quality. If you can't solve the abnormal situation by yourself, you can use electricity</p><p>Consult the manufacturer or the unit designated by the manufacturer.</p></div><div><div>Enclosure</div><div><p>1 manual</p></div></div></div>
---	---	---